

FOR DEMONSTRATION PURPOSES ONLY - NOT FOR CLINICAL USE

# SAMPLE COLLECTION INSTRUCTIONS

## YOUR KIT INCLUDES




## STEP 1: ADULT PREPARES THE KIT

**1.1** Wash or sanitize your hands and **PLACE ALL THE CONTENTS** of the kit on a clean dry surface.




Ensure the child has washed their hands.

**1.2** Use a permanent marker to **FILL IN TWO ID STICKERS WITH THE REQUESTED DATA.** Then stick one completed ID sticker onto the tube.

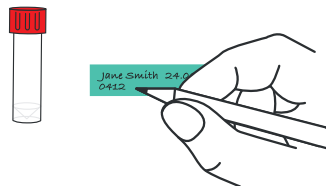


**1.3 TOUCHING ONLY THE HANDLE,** remove the swab out of its packaging.

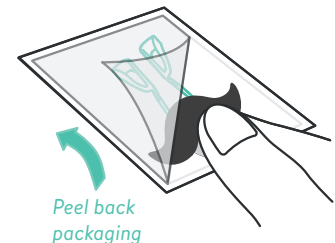


## STEP 2: CHILD SELF-COLLECTS THE SAMPLE

**2.1 ASK THE CHILD** to hold the swab by the handle and bring it to their nose.

**STEP 2** should be supervised by an adult.


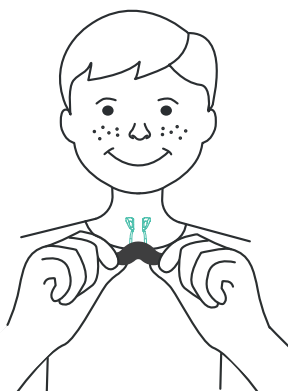


**2.2 THE CHILD SHOULD PLACE THE SWAB** gently in their nose as far as is comfortable.

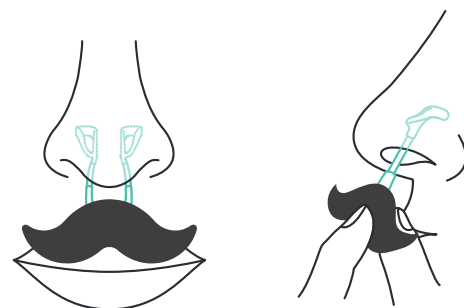

### 2.3 SET TIMER FOR 15 SECONDS

and ask the child to breathe normally through their nose.

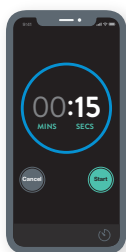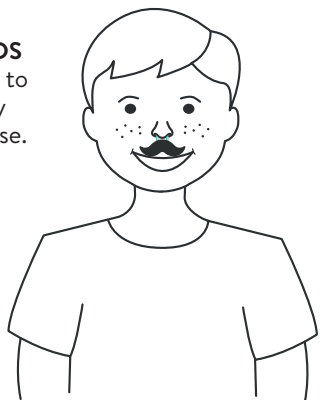

Set timer  
for 15 SECS


**2.4 AFTER 15 SECONDS,**  
ask the child to move the swab  
inside their nose for **15 seconds**  
making sure the loops make  
continual contact with the  
inside of their nose.



## STEP 3: ADULT PREPARES THE SAMPLE

### 3.1 BREAK BOTH LOOPS INTO THE TEST TUBE ONE AT A TIME.


Holding the tube firmly in one hand and the swab in the other, insert **ONE LOOP AT A TIME** into the test tube and snap off at the break line.

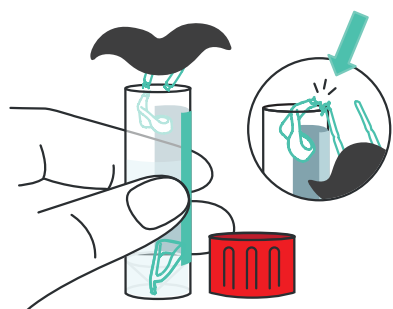

Do this by bending the swab arm against the edge of the tube at the break line.

### 3.2 SCREW THE CAP BACK TIGHTLY ON THE TUBE

ensuring it is tightly sealed so the liquid does not leak.

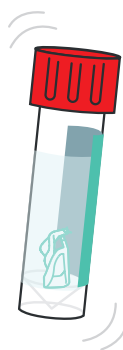


SHAKE THE  
TUBE GENTLY A  
FEW TIMES

Put the swab handle  
in the refuse bag  
and discard.




### 3.3 USE AN ALCOHOL WIPE

to sanitize around the top  
of the tube.

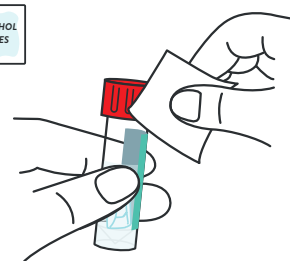

### 3.4 PLACE THE TUBE INTO THE BIOHAZARD BAG and sanitize the bag.


Discard the alcohol wipe.

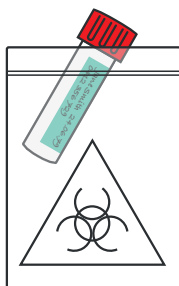

### 3.5 CLEAN YOUR HANDS

Apply hand sanitizer, wash your hands thoroughly or use the alcohol wipe provided.

Discard the alcohol wipe.

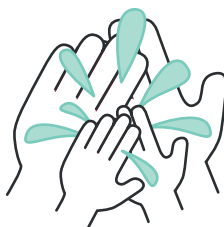

## STEP 4: RETURN THE SAMPLE

Please follow the instructions to return the sample to the relevant testing center.



HOSPITAL • FEVER CLINIC • LAB

